

# Galaxy Fixation

Gemini™ System



# Galaxy Fixation

Gemini™ System

## Double Design, Infinite Solutions

*The Galaxy Fixation Gemini™ System is a modular external fixation system for fracture treatments of lower and upper limbs. It can be used as a hybrid system in conjunction with Orthofix circular external fixators.*

### Simple

Single clamps can perform rod-to-rod coupling and pin-to-rod coupling.

### Versatile

Universal clamps can accommodate different diameters of rods and bone screws, providing the modular system with mounting versatility for each anatomical district.

### Stable

Internal teeth and locking profiles designed to provide high torsional strength and avoid components sliding.

### Easy to use

Snap-in system, provisional tightening by hand, definitive cam closure in one step. On each clamp slot there is a marker to identify proper rod/bone screw diameter. Preassembled multipin clamps.



# Galaxy Fixation Gemini Universal Clamps

In a multiple use non-sterile configuration

Markers for easy identification of rod and pin diameter

Modified knurled knob and definitive closure mechanism

Intuitive design with visible check for definitive locking regardless of clamp position

Hand-tightening pre-closure

Notched interface profile of the clamp and wrench definitive tightening

- Simplify intra-operative reduction adjustment
- Fast frame assembly
- Ease of use because there is no need for instruments for primary frame stabilization/closure

Provide adequate stability of the frame after final closure


Available also in multipin version

**Galaxy Fixation Gemini  
UNIVERSAL SINGLE CLAMP  
(94100)**

**Galaxy Fixation Gemini  
DOUBLE UNIVERSAL MULTIPIN  
CLAMP (94200)**

**Galaxy Fixation Gemini  
UNIVERSAL MULTIPIN  
CLAMP (94300)**

Allows to perform 10 different types of coupling

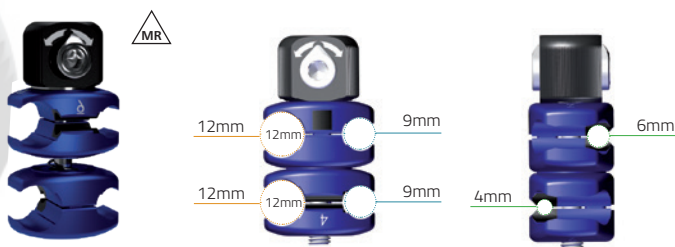





# Galaxy Fixation Gemini Standard Clamps

In single use sterile configuration

Galaxy Fixation Gemini  
SINGLE CLAMP  
(99-94010)

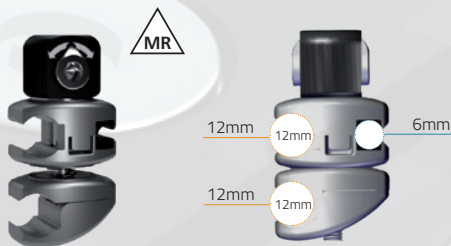

Galaxy Fixation Gemini  
TRANSITIONAL SINGLE  
CLAMP (99-94030)

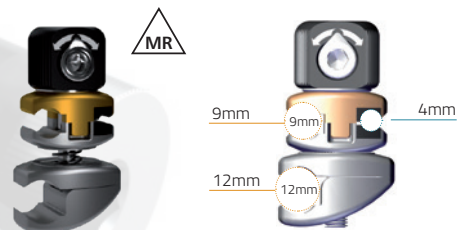

Available also in multipin version

Galaxy Fixation Gemini  
DOUBLE MULTIPIN CLAMP  
MEDIUM (99-94140)



Galaxy Fixation Gemini  
DOUBLE MULTIPIN CLAMP  
(99-94040)



Galaxy Fixation Gemini  
MULTIPIN CLAMP  
(99-94020)






All these clamps are provided sterile.

# Galaxy Fixation Gemini™ System

## All-In-One Sterilization Tray

Galaxy Fixation Gemini™ Universal Clamps are multiple use and available in non-sterile configuration

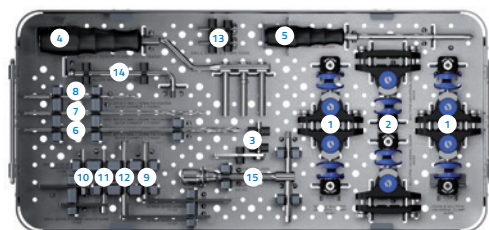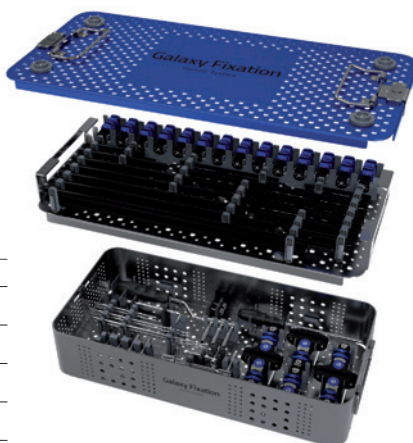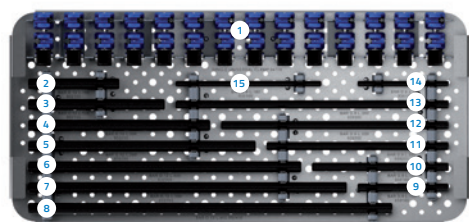

Part #	Description	Quantity
1 94100	GALAXY FIXATION GEMINI UNIVERSAL SINGLE CLAMP	14
2 932100	ROD D 12 MM L 100 MM	2
3 932150	ROD D 12 MM L 150 MM	2
4 932200	ROD D 12 MM L 200 MM	2
5 932250	ROD D 12 MM L 250 MM	2
6 932300	ROD D 12 MM L 300 MM	2
7 932350	ROD D 12 MM L 350 MM	2
8 932400	ROD D 12 MM L 400 MM	2
9 939100	ROD D 9 MM L 100 MM	2
10 939150	ROD D 9 MM L 150 MM	2
11 939200	ROD D 9 MM L 200 MM	2
12 939250	ROD D 9 MM L 250 MM	2
13 939300	ROD D 9 MM L 300 MM	2
14 936100	ROD D 6 MM L 100 MM	2
15 936160	ROD D 6 MM L 160 MM	2

Part #	Description	Quantity
1 94200	GALAXY FIXATION GEMINI DOUBLE UNIVERSAL MULTIPIN CLAMP	2
2 94300	GALAXY FIXATION GEMINI UNIVERSAL MULTIPIN CLAMP	2
3 19995	SCREW GUIDE	2
4 194500	GALAXY FIXATION GEMINI SCREW GUIDE	1
5 194400	GALAXY FIXATION GEMINI TROCAR	1
6 1-1100201	DRILL BIT D.4.8 MM L.240 MM TIN COATED - QUICK CONNECT	2
7 1-1100301	DRILL BIT D.3.2 MM L.200 MM TIN COATED - QUICK CONNECT	2
8 1-1355001	DRILL BIT D.2.7 MM L.127 MM TIN COATED - QUICK CONNECT	2
9 11102	SCREW GUIDE, LENGTH 60 MM	2
10 11106	DRILL GUIDE DIAMETER 3.2 MM LENGTH 40 MM	2
11 11137	SCREW GUIDE, LENGTH 80 MM	2
12 11138	DRILL GUIDE DIAMETER 4.8 MM LENGTH 60 MM	2
13 13530	DRILL GUIDE DIAMETER 2.7 MM	2
14 30017	ALLEN WRENCH 5 MM	1
93162 or 91150	T-WRENCH HEXAGON 5-5 QC or UNIVERSAL T-WRENCH	2

### Galaxy Fixation Gemini UNIVERSAL CLAMPS (94100, 94300)

are fully compatible and can be used with

### Galaxy Fixation ELBOW HINGE (93410)




# Galaxy Fixation Gemini™ System


## Ready to Use Sterile Sets

Galaxy Fixation Gemini™ Standard Clamps are single use and available in sterile sets or single sterile packs

Please refer to GF-2101-OPT Galaxy Fixation Gemini™ Operative Technique for complete overview

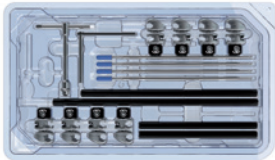


99-94710 Gemini Pelvis Sterile Set

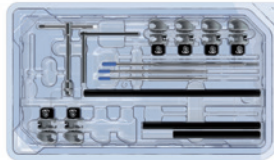


99-94730 Gemini Ankle 6mm Sterile Set



99-94750 Gemini Metatarsal Sterile Set\*

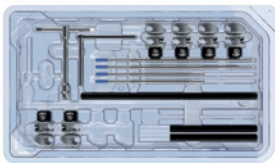


99-94720 Gemini Z Frame Sterile Set

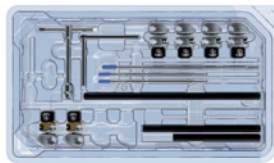


99-94740 Gemini Ankle 4mm Sterile Set



99-94004 Gemini Foot Support Sterile Set\*

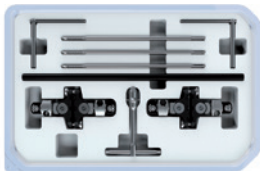

99-94760 Gemini Quadrilateral Sterile Set

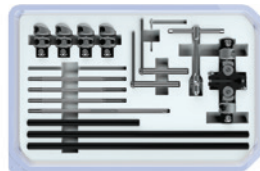

99-94001 Gemini Ankle 6mm Multiscrew Clamp Sterile Set

\*To be used in combination with Ankle Set (99-94730 or 99-94740 or 99-94001)

Please refer to the "Instructions for Use" supplied with the product for specific information on indications for use, contraindications, warnings, precautions, adverse reactions and sterilization.

Electronic Instructions for use available at the website <http://ifu.orthofix.it>

Electronic Instructions for use - Minimum requirements for consultation:

- Internet connection (56 Kbit/s)
- Device capable to visualize PDF (ISO/IEC 32000-1) files
- Disk space: 50 Mbytes

Free paper copy can be requested from customer service (delivery within 7 days):

tel +39 045 6719301, fax +39 045 6719370,

e-mail: [customerservice@orthofix.it](mailto:customerservice@orthofix.it)

Caution: Federal law (USA) restricts this device to sale by or on the order of a physician. Proper surgical procedure is the responsibility of the medical professional. Operative techniques are furnished as an informative guideline. Each surgeon must evaluate the appropriateness of a technique based on his or her personal medical credentials and experience.



Manufactured by:  
ORTHOFIX Srl  
Via Delle Nazioni 9, 37012 Bussolengo  
(Verona), Italy  
Telephone +39 045 6719000, Fax +39 045  
6719380

CE 0123

**Rx Only**


Orthofix products or services referenced herein are trademarks or registered trademarks of Orthofix SRL and its group of companies. All rights reserved.

**Orthofix.com**  
**GalaxyFixation.com**



GF-2105-PL-E0 C 05/23

Distributed by:

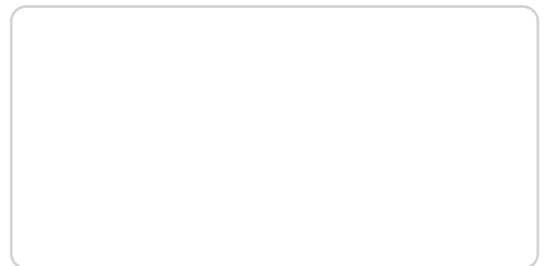

 **ORTHOFIX**<sup>®</sup>  
ORTHOPEDICS