

TrueLok™ Elevate

Transverse Bone Transport System
Sales Sheet



The Strength of Tradition,
the Innovation of Regeneration



Built on a proven foundation, designed to treat bony and soft tissue defects and deformities with state-of-the-art clinical approach

The TrueLok™ Elevate is a sterile-packed external fixator indicated for the treatment of:

- Bony and soft tissue defects
- Bony and soft tissue deformities
- Non-unions (U.S. Only)

It allows for transverse distraction and compression movement of the bone segment.


It can be used as a stand-alone device or in conjunction with TrueLok™, TL-HEX™ and TrueLok™ EVO.

99-50-2513 TRUELOK ELEVATE SET STERILE

Part Number	Description	Q.ty
1-1355001	Drill Bit D.2,7mm L.127mm TIN	1
50-25071	Drill Insert	2
50-25070	Wire Insert	2
50-2505	TrueLok Elevate External Fixator	1
50-2506	TrueLok Elevate Template	1
50-2511	Half Pin L120mm D4.0mm Thread D4.0xL18mm	2
50-2508	K-Wire L180mm D1.8mm	3


50-2590C TRUELOK ELEVATE INSTRUMENT TRAY COMPLETE

Part Number	Description	Q.ty
M210	T-Wrench For Bone Screws	1
50-2510	TL Straight Tube Wrench	1
54-1154	TL, Wrench, Combo, 10mm	1
93162	T-Wrench Hexagon 5-5 QC	1
54-2226	TL, 90 Degree Tubular Wrench	1
50-2509	TrueLok Elevate Drill Positioning Guide	1
W1003	Wire Cutter	1
54-2227	Needle Nose Pliers, Stainless Steel	1
11104	Drill Guide 4.8mm L40mm	1
50-2515	Drill Guide 3.2mm L40mm	1

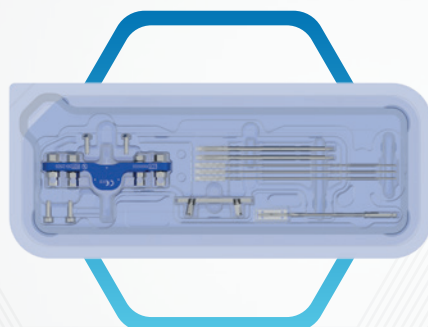


Simple

The fixator is pre-assembled: it allows a quick fixation of the device with the half-pins

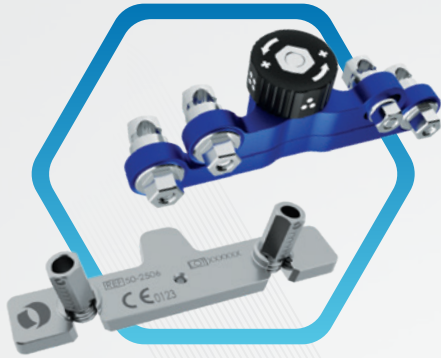


Versatile

Design allows device to be attached even if half pins are in different planes

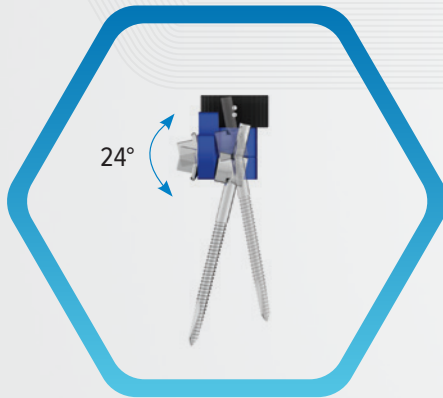

Stable

Half pin fixation bolts with serrations and grooves provide superior bone segment stabilization

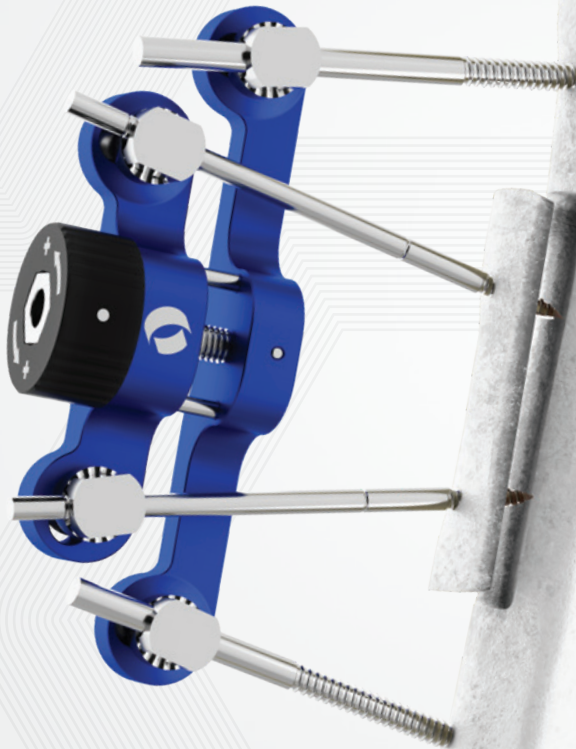


Reproducible and Standardized Operative Technique

One Mechanism to Activate the Device

The specifically designed knob allows for 0.25mm lateral distraction-compression each quarter of turn

0.25mm/quarter of turn

Running distance up to 14mm

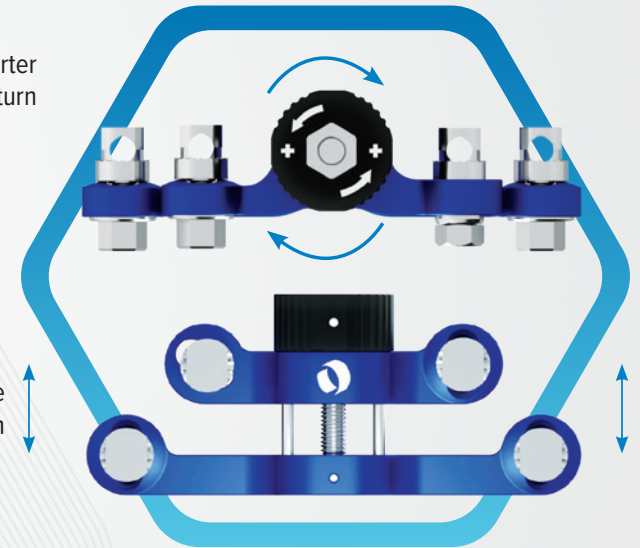

Time Saving

All the components are provided sterile: no maintenance/sterilization cost for the hospitals


**Manufactured by:**

ORTHOFIX Srl
Via Delle Nazioni 9,
37012 Bussolengo (Verona), Italy
Tel.: +39 045 6719000
Fax: +39 045 6719380

CE₀₁₂₃ **Rx Only**

Orthofix products or services referenced herein are trademarks or registered trademarks of Orthofix SRL and its group of companies. All rights reserved.

TrueLok™, TL-HEX™, and TrueLok™ EVO are trademarks or registered trademarks exclusively licensed to Orthofix Medical Inc., and its affiliate companies. All rights reserved.

Distributed by:

orthofix.com

TL-2407-PL-E0 B 01/25

